

# Zimmer® Dermatome AN

## Meeting Your Skin Grafting Needs

With more than 50 years of experience in skin graft instrument design and manufacturing, Zimmer Biomet offers the Zimmer Dermatome AN. This instrument provides better results with consistent, reliable skin graft retrieval.

Additional features include:

- Autowashable in neutral pH cleaners/detergents
- Safe handling with proprietary blade and blade cover
- Sealed to enable soaking and easy removal of debris
- Width plate secured in place in three locations
- Easy assembly with one width plate screw
- Cutting range from 0.004" - 0.024" (0.1 mm – 0.60 mm)



## Easy Maneuverability. Easy Maintenance.

Experience power and durability without compromising maneuverability. This lightweight device is easy to hold and control with less fatigue for the surgeon. Designed for easy maintenance, the device system does not require lubrication.

# Feature Details

## Handpiece Features

### Width Plate Tabs

- Secures the width plate to the dermatome
- Ensures proper assembly
- Holds width plate in place during retrieval



### One Width Plate Screw

- Enables quicker set up/blade change
- Provides efficient use of OR time



### Laser Etching

- Easier to read adjustment knob



### Sealed Unit

- Protects motor from water ingress and corrosion
- Enables soaking for dried on debris
- Allows auto washing in neutral pH detergent



### Anodized Coating

- Allows for auto washing device
- Neutral pH detergent (6.6 – 7.3 on a pH scale)

## Service and Repair Programs

All Zimmer Biomet dermatomes are recommended for service and electronic calibration annually. Learn more at [www.service.zimmer.com](http://www.service.zimmer.com).

An annual Clinical Excellence visit is included with Premium Warranty on the Zimmer Dermatome AN during the duration of the service agreement. For more information contact [SurgicalSales-ClinicalExcellence@zimmerbiomet.com](mailto:SurgicalSales-ClinicalExcellence@zimmerbiomet.com)

This documentation is intended exclusively for physicians and is not intended for laypersons. Information on the products and procedures contained in this document is of a general nature and does not represent and does not constitute medical advice or recommendations. Because this information does not purport to constitute any diagnostic or therapeutic statement with regard to any individual medical case, each patient must be examined and advised individually, and this document does not replace the need for such examination and/or advice in whole or in part. Please refer to the package inserts for important product information, including, but not limited to, indications, contraindications, warning, precautions, and adverse effects.


©2017 Zimmer Biomet



## Blade Features

### Blade cover

- Easy to remove after blade is properly inserted
- Minimizes sharps exposure when inserting/removing the blade

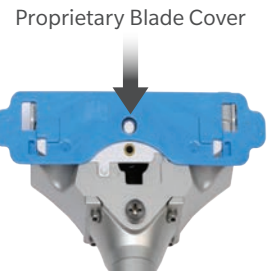

### Notch and Key

- Mates with head design for proper assembly
- Lowers the chance of incorrect assembly
- Reduces the risk of patient injury due to improper assembly



### Brass Ring

- Prevents drive pin from enlarging the hole in the blade
- Ensures full travel distance of the blade for optimal cutting



### Finger Cut Outs

- Allows easy access to lift off blade
- Improves staff safety



Contact your Zimmer Biomet representative or visit us at [www.zimmerbiomet.com/surgical](http://www.zimmerbiomet.com/surgical)

## Ordering Information

Description	Part No.
Zimmer Dermatome AN	88710100
Zimmer Dermatome Hose	00880100200
Autoclave Case	88700300

### Accessories

Z Dermatome Blades	88700010
Zimmer Dermatome AN Screwdriver	88700500
1 in. (2.5 cm) Width Plate	8870N100
1.5 in. (3.8 cm) Width Plate	8870N150
2 in. (5.1 cm) Width Plate	8870N200
3 in. (7.6 cm) Width Plate	8870N300
4 in. (10.2 cm) Width Plate	8870N400
Width Plate Screws, 5 PK	88701805
Zimmer Dermatome AN Manual	88711FU
O-Rings	00880110541