


The image displays two teal-colored metal rods of the FastFrame™ External Fixation System. Each rod is equipped with black clamps at both ends, which are used to secure the rods to a patient's bone. Red wing-like handles are attached to the clamps for manual adjustment. The rods are positioned diagonally across the frame. A blue banner is overlaid on the left side, and a light blue banner is overlaid in the center. The Zimmer Biomet logo is visible on the lower part of the rods.

FastFrame™ External Fixation System

Product Brochure

 ZIMMER BIOMET
Your progress. Our promise.®

High Quality, Low Cost Materials

FastFrame kits are pre-assembled, sterile packed, and 100% disposable; meeting the demands of both surgeon and hospital.

- Simplicity and rapidity are the design philosophies
- May eliminate risk of tray contamination
- Reduces need for reprocessing and re-sterilization
- Reduce logistical burden to the hospital
- High quality, pre-assembled materials maintain frame integrity





Product Features

The FastFrame External Fixation System removes the complexity, lengthiness, and potentially excessive costs from common external fixation procedures.

These new external fixation frames introduce several new advancements:




Telescoping Tubes

Telescoping Tubes (1) eliminate the need for multiple bar lengths and are bi-directional allowing independent distraction and compression. Pull the red tabs (2), and the frame becomes uni-directional allowing for easy distraction and resists compression.

Tri-Planar Alignment Spheres

The spheres allow for almost limitless pin placement, alignment, and rotation (3). The rigid frame is designed for efficient and quick application. Each frame provisionally locks in place via wing nuts (4).



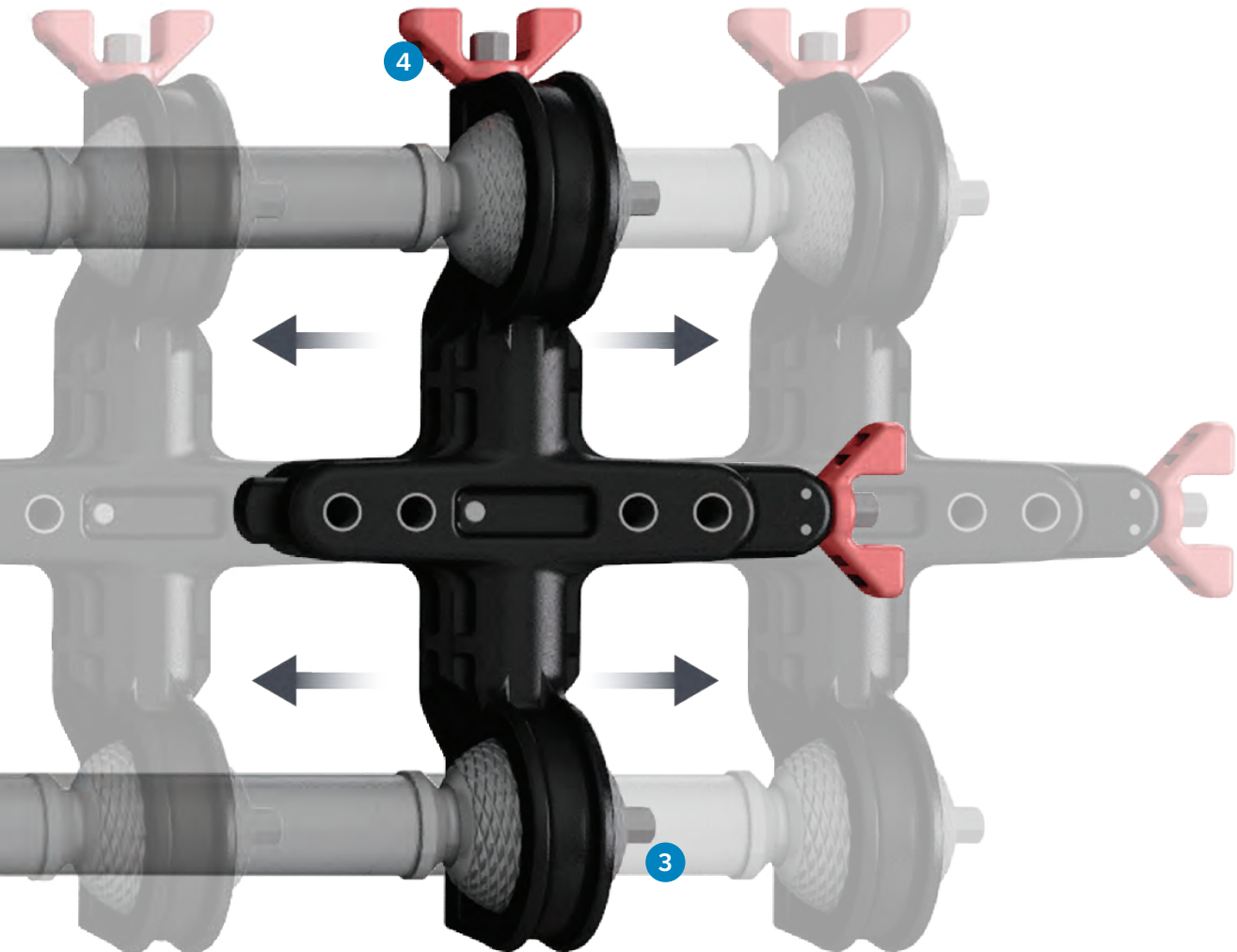

Applications

FastFrame sterile kits are being offered for the ankle spanning, distal radius, damage control, and knee spanning applications.



Ankle Spanning



Distal Radius



Damage Control



Knee Spanning

All content herein is protected by copyright, trademarks and other intellectual property rights, as applicable, owned by or licensed to Zimmer Biomet or its affiliates unless otherwise indicated, and must not be redistributed, duplicated or disclosed, in whole or in part, without the express written consent of Zimmer Biomet.

This material is intended for health care professionals. Distribution to any other recipient is prohibited.

For product information, including indications, contraindications, warnings, precautions, potential adverse effects and patient counseling information, see the package insert and www.zimmerbiomet.com.

Check for country product clearances and reference product specific instructions for use.

Not for distribution in France.

©2018 Zimmer Biomet



ZIMMER BIOMET
Your progress. Our promise.®

0158.1-GLBL-en-REV0318



Zimmer, U.K. Ltd.
9 Lancaster Place
South Marston Park
Swindon, SN3 4FP, UK

www.zimmerbiomet.com

